


LANE 1 2 3 4 5 6 7 8

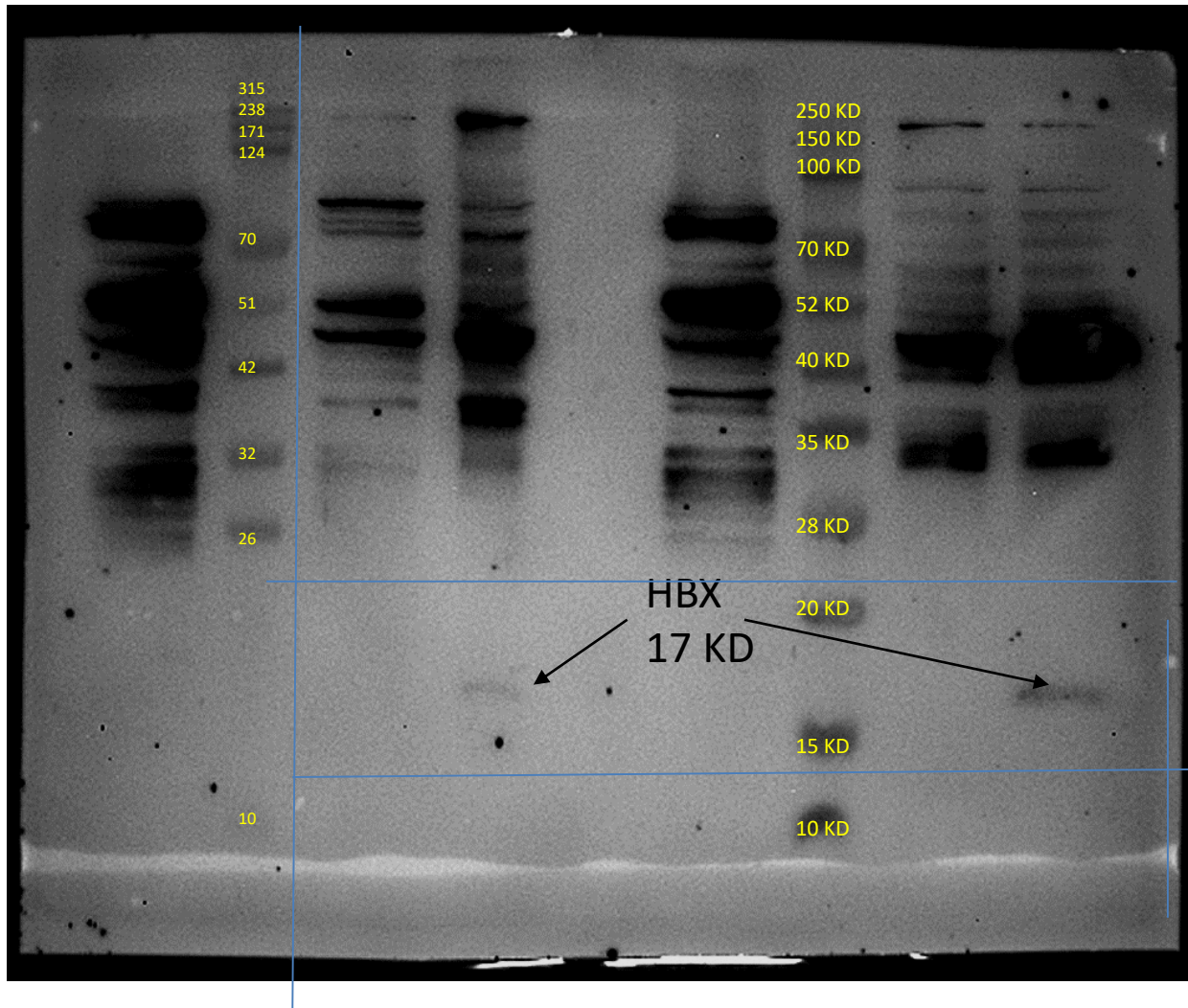

LANE

1

2

3

4

5

6

7

8

SAMPLE

ES

MARKER 1

Fibroblast control


Fibroblast Tx pHBx

ES

HEK Control

HEK Tx pHBx

Figure-1A

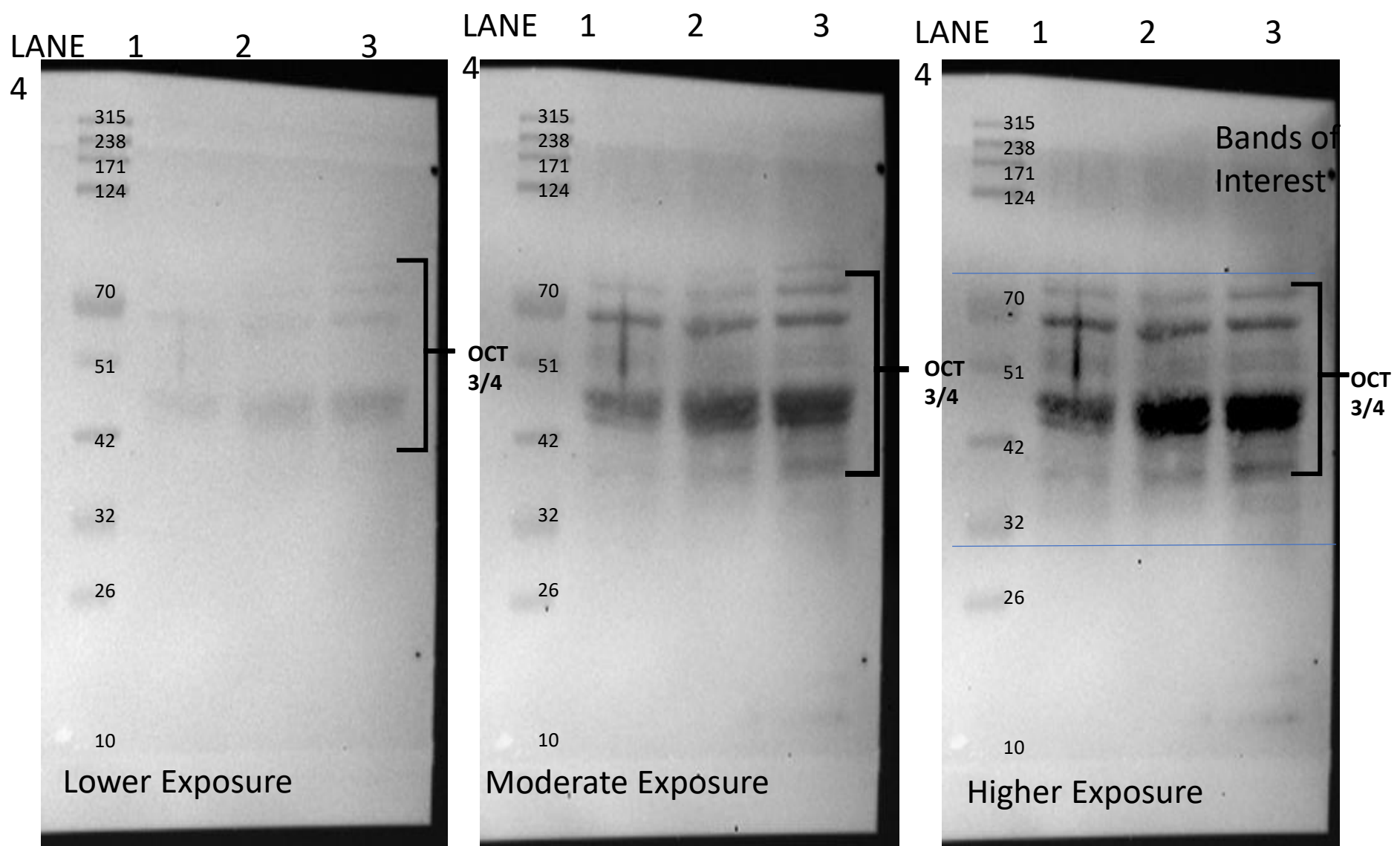

Lane	Sample
1	IPSC(10 ul)
2	IPSC(15 ul)
3	IPSC(20 ul)


MW of OCT-3/4A 52 KD
 MW of OCT-3/4B 45 KD

I have marked the total band position referring to the datasheet of the antibody(SC-5279). Although I have marked separately our bands of interest as far as my understanding considering the multiple isoforms of OCT based on available literature.

Figure-D-g

Lane 0 1 2 3 4 5 6 7 8

Figure-1-E-a
W/B

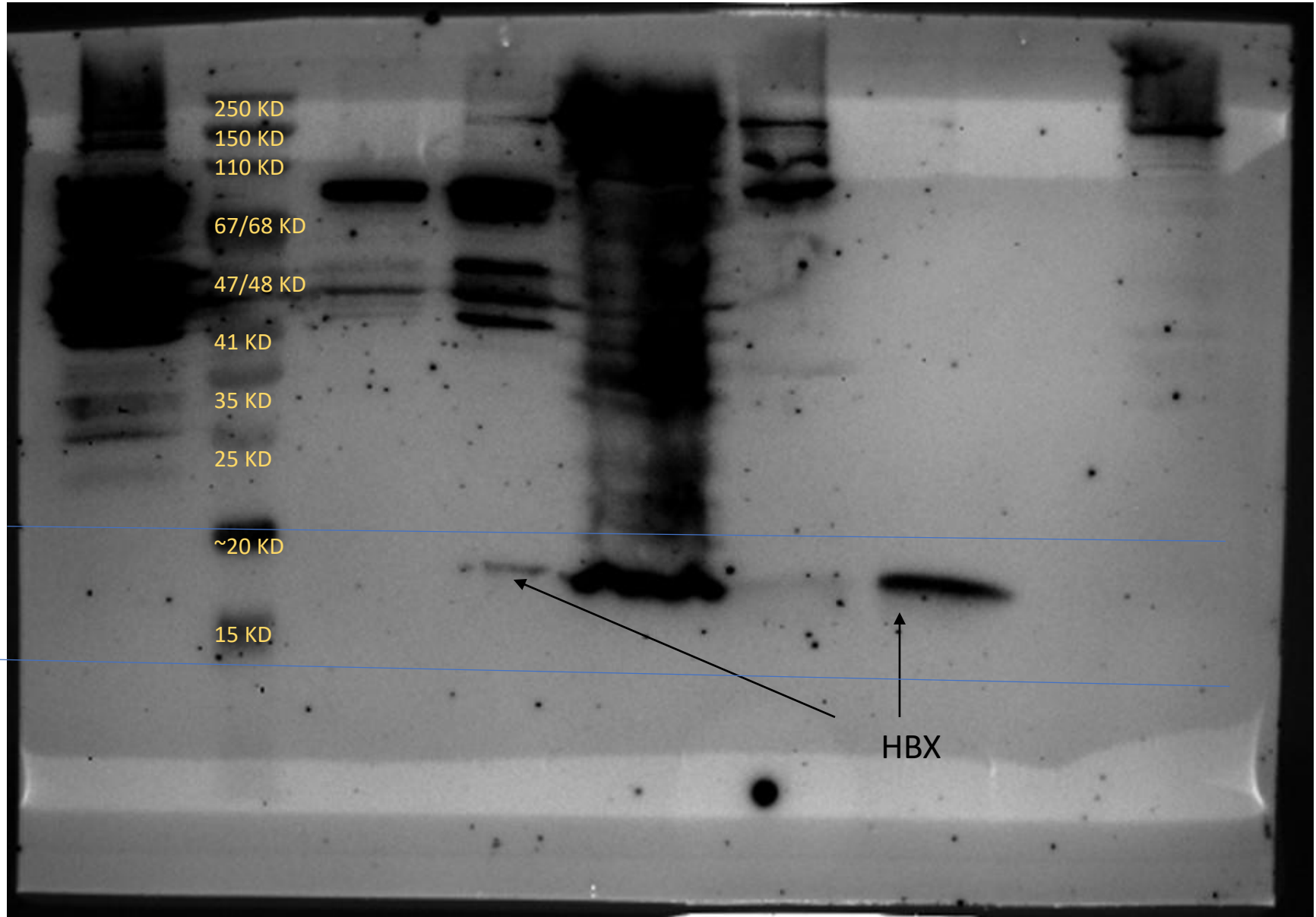

Lane	Sample
0	HBx-iPSC
1	Marker
2	Fib control
3	Fib pHBX TF
4	HEK pHBX TF 5mcg
5	NIIHEKHBX TF 0.5mcg
6	HEK HBX TF 2 mcg
7	HEK Control
8	HEK Control-3x

12% GEL
Blocked with 5% BSA
Primary antibody:
Hepatitis B Virus X
Monoclonal Antibody
(X36C)[MA1-081]: 1:1000
Secondary ab: anti-
mouse:1:10000

Data for HBX expression
in HEK and Fibroblast


Marker pHbX iPSC

Figure-1-E-b
(DNA gel)



0.5kb

Lane 1 2 3 4 5

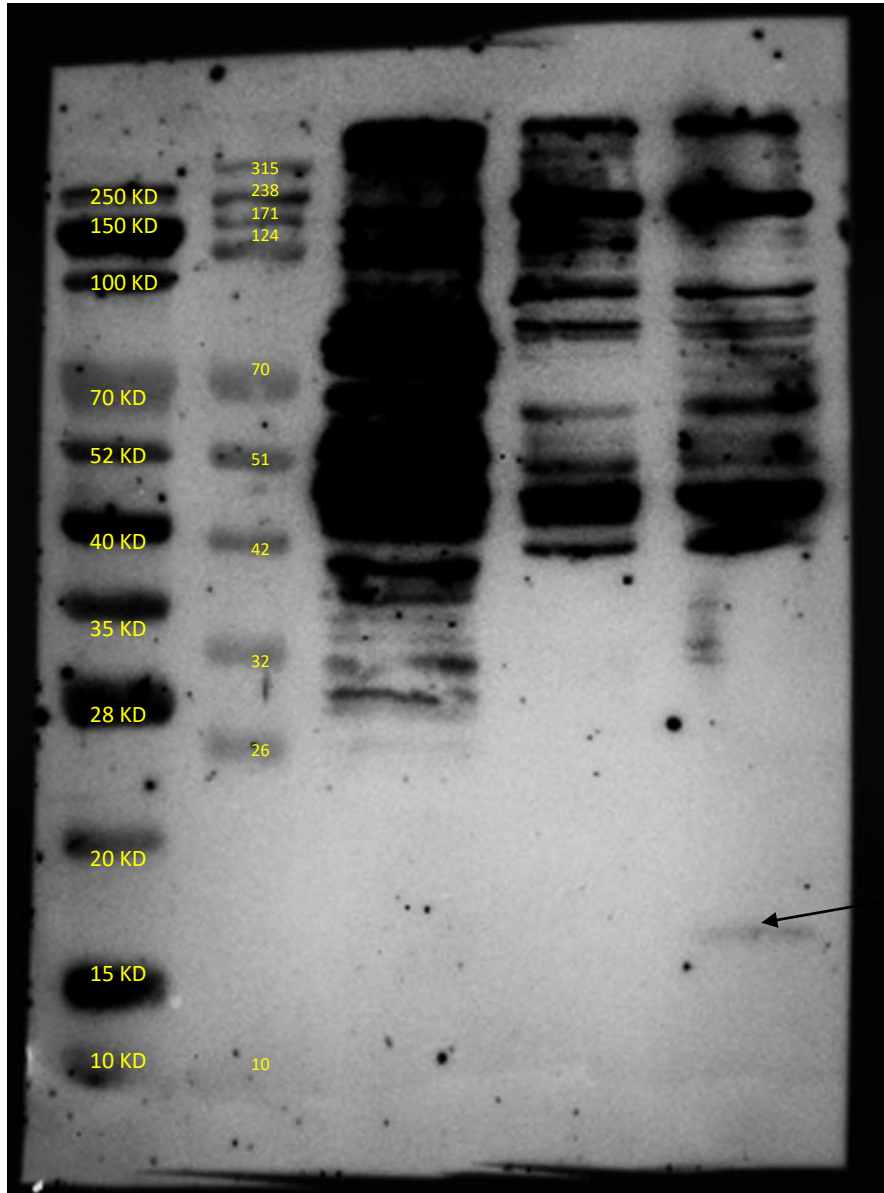


Lane	Sample
1	Marker 1
2	Marker 2
3	IPSC
4	HEK Control
5	HEK pHBX TF


HBX
17KD

Figure-5A
Row-1

Lane 1 2 3 4 5

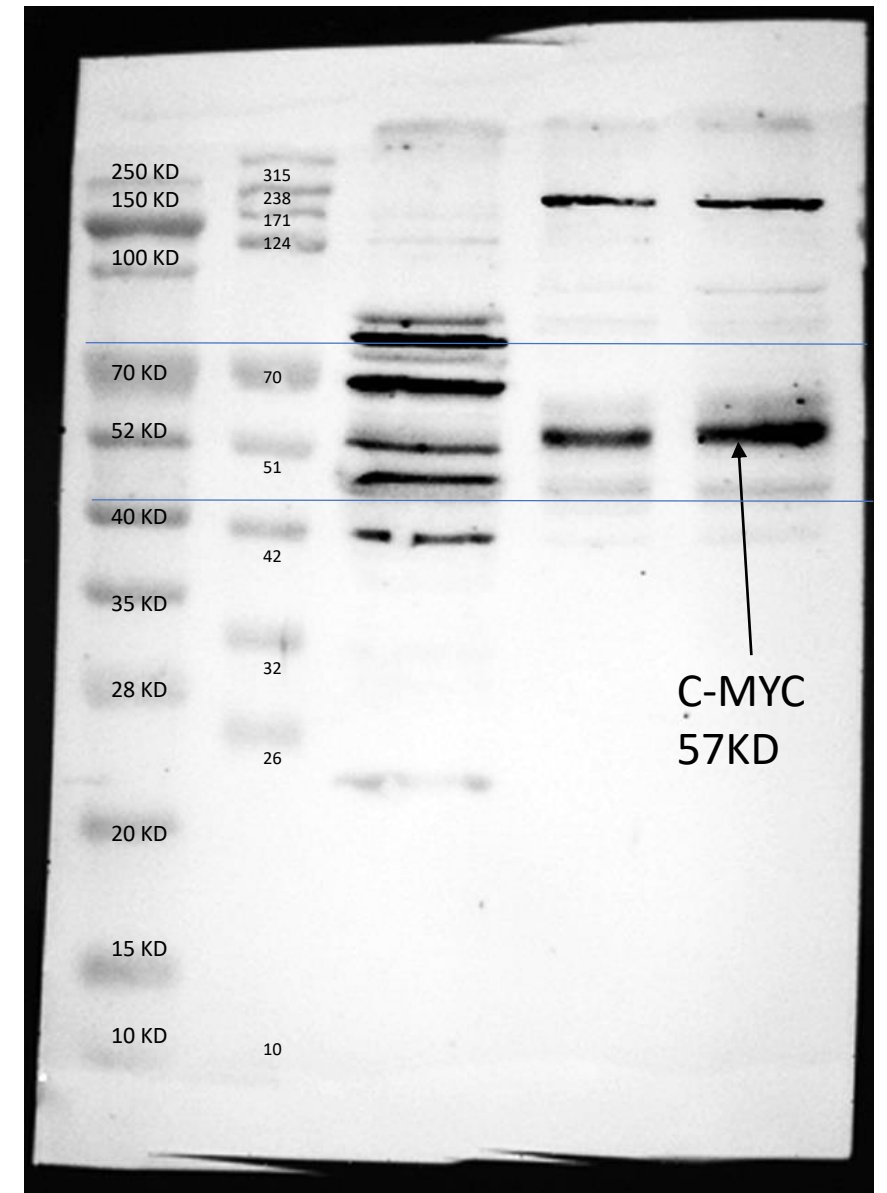


C-MYC
57KD

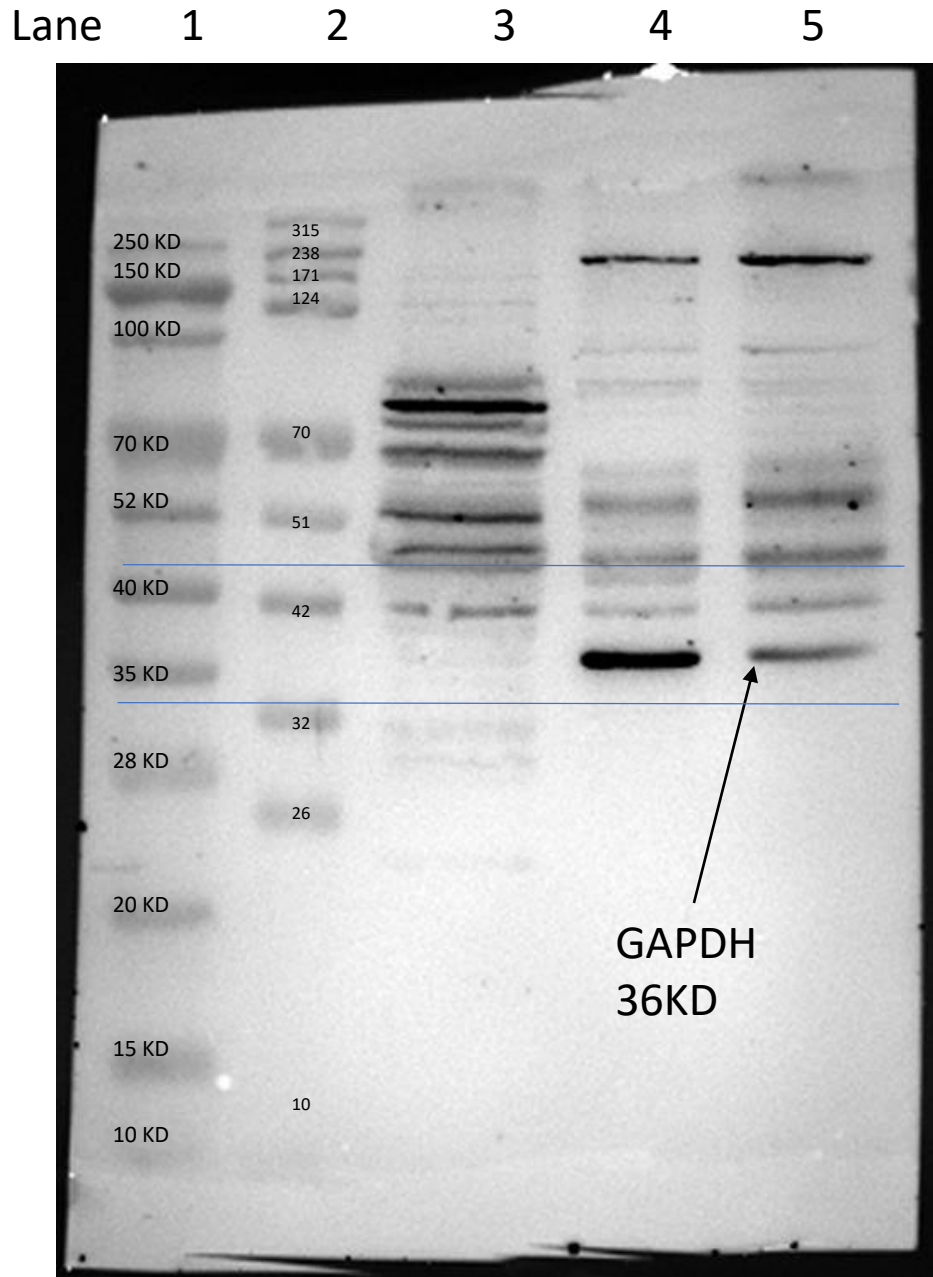

Lane	Sample
1	Marker 1
2	Marker 2
3	IPSC
4	HEK Control
5	HEK pHBX TF

HBX
↓
C-MYC
↓
GAPDH

Figure 5A
Row 2

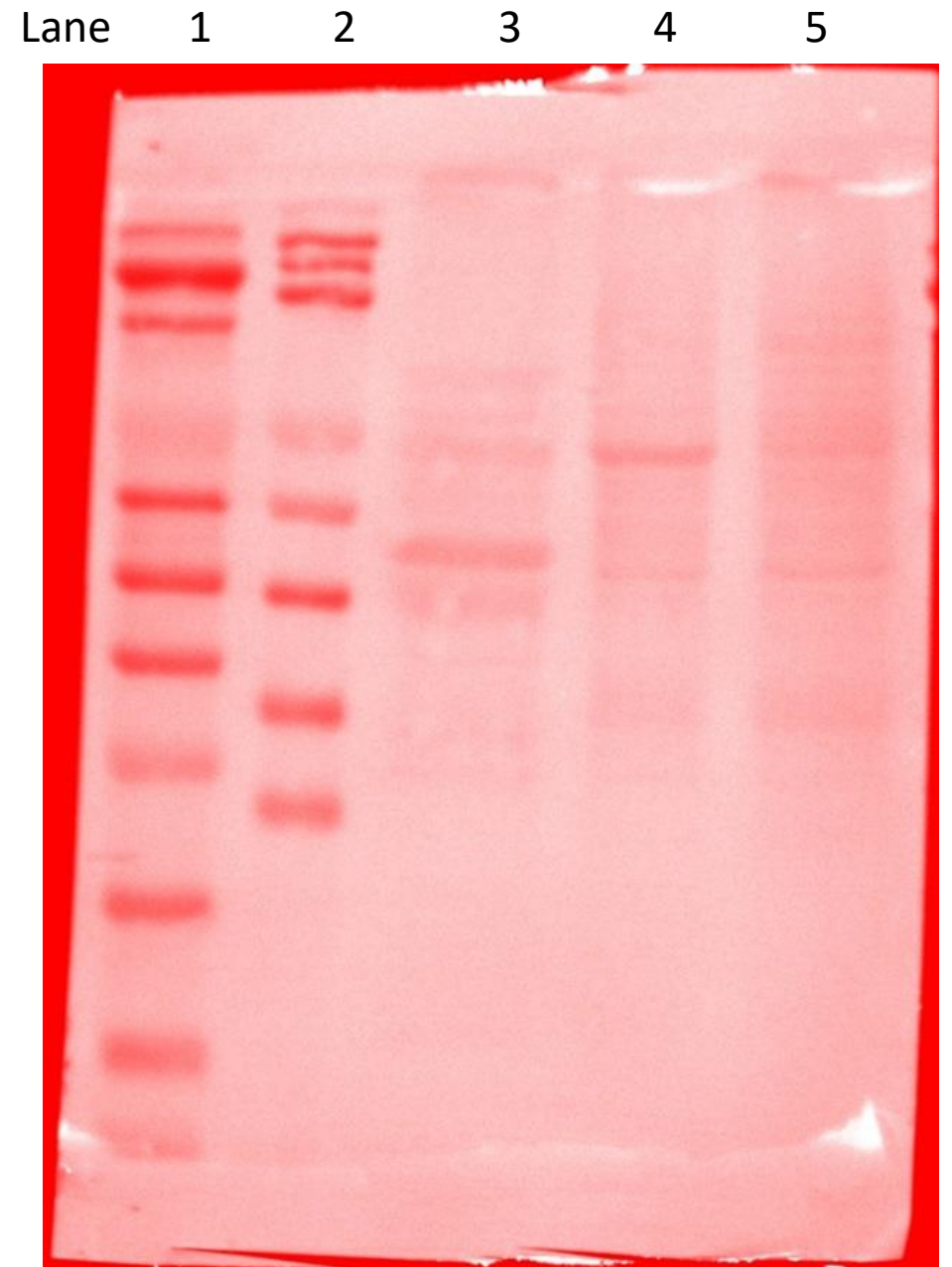





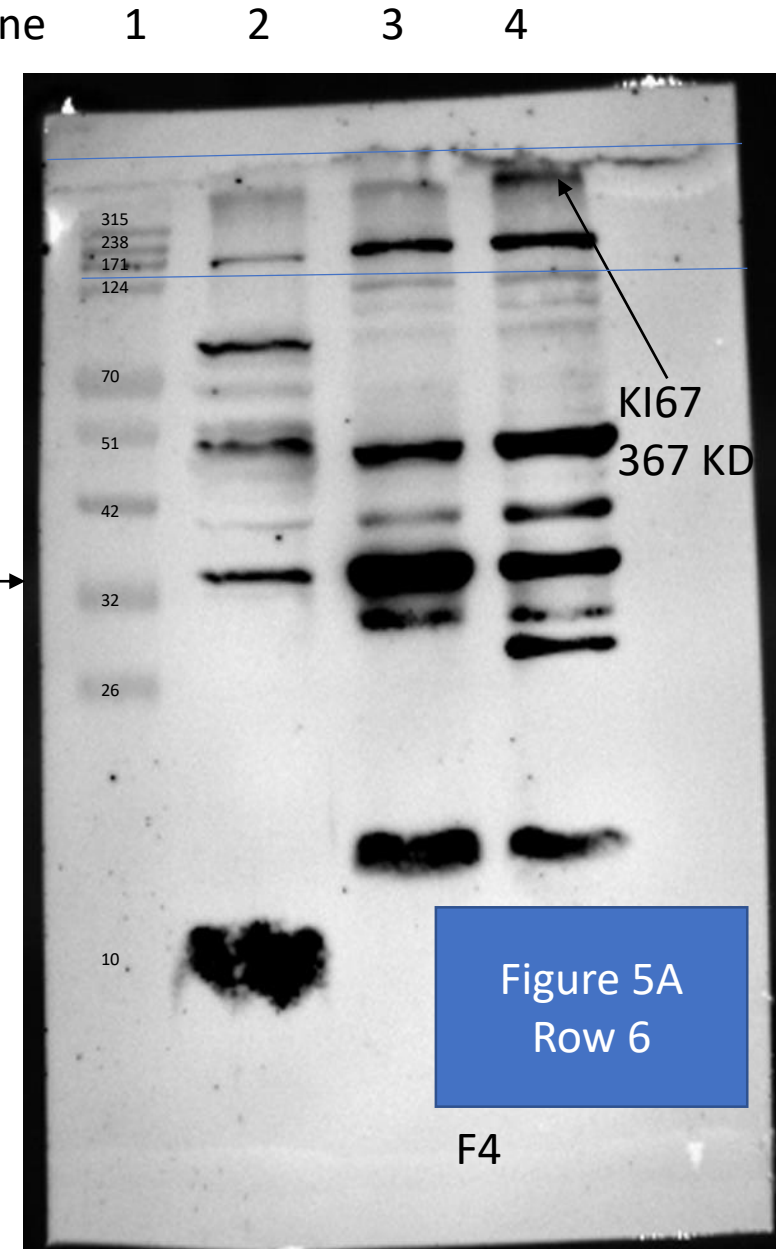


Figure 5A
Row 3

LANE SAMPLE
1 IPSC
2 HEK CONTROL
3 HEK pHBX TF

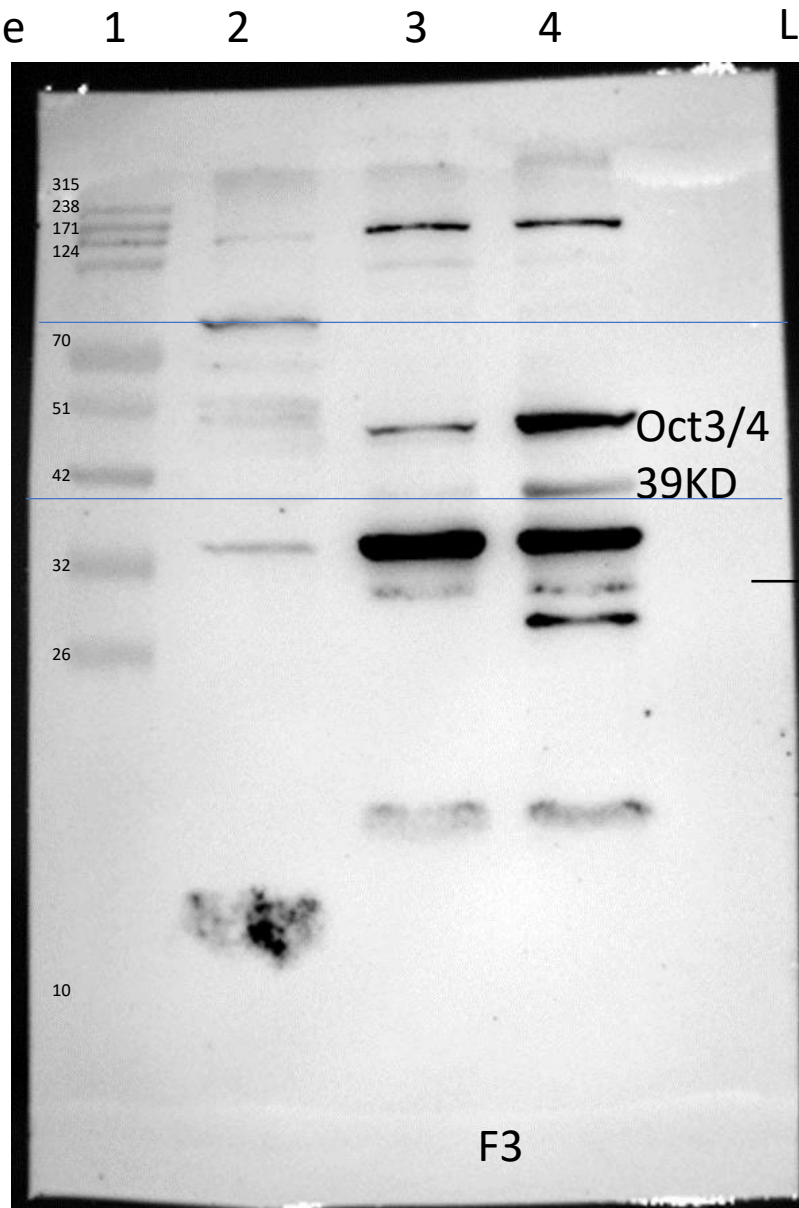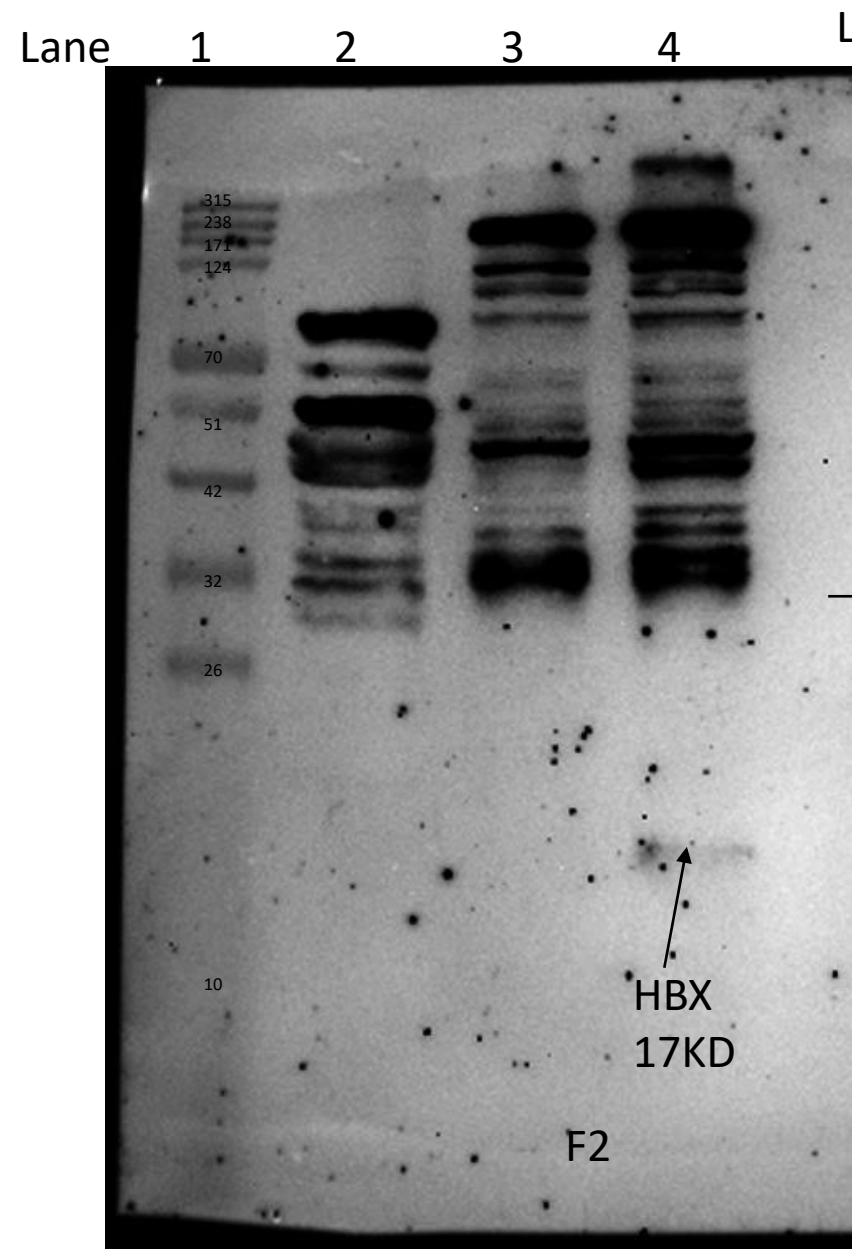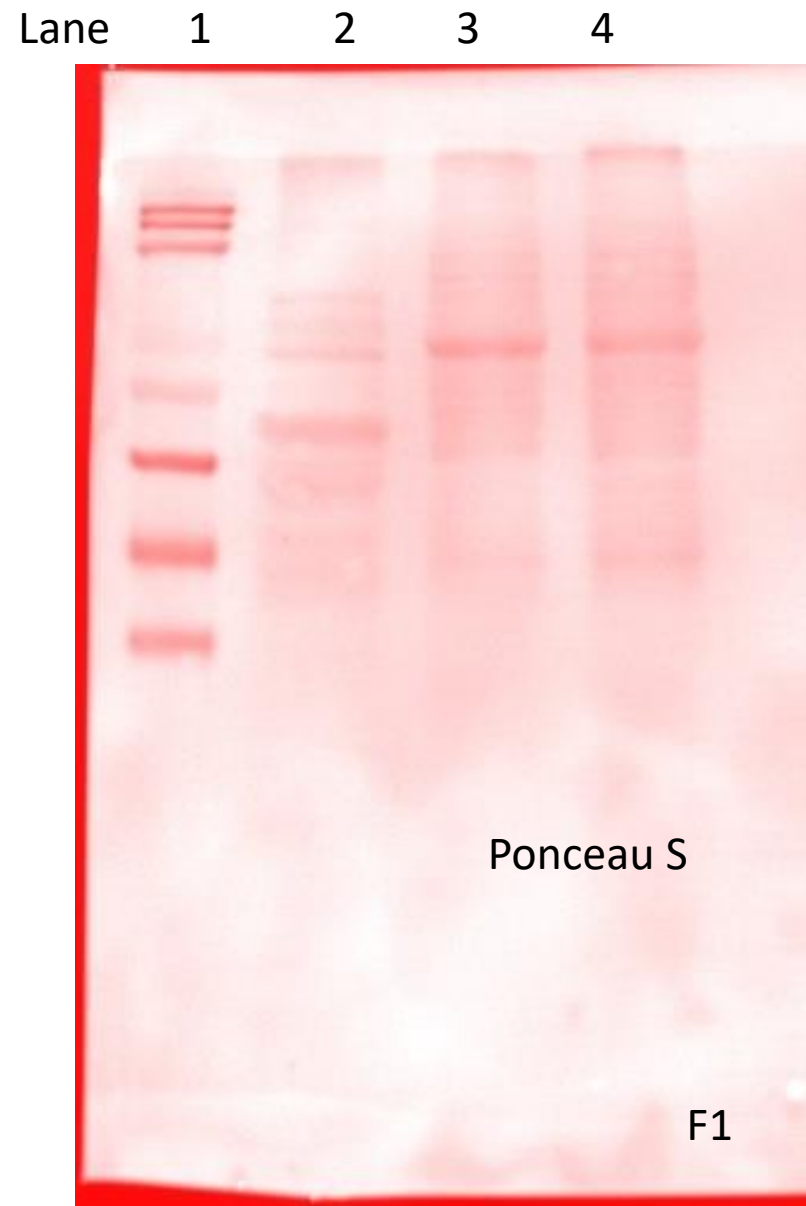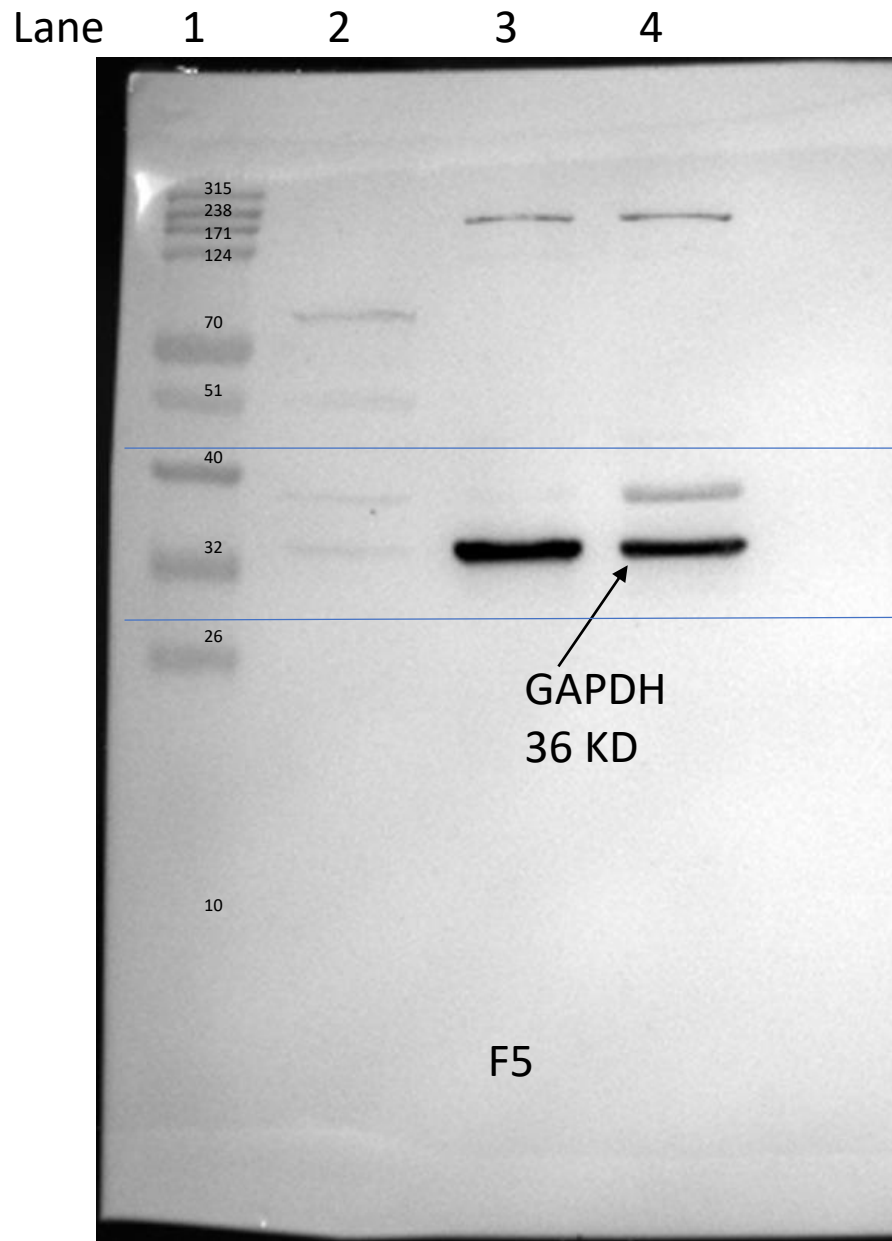


LANE	SAMPLE
1	IPSC
2	HEK CONTROL
3	HEK HBX TF

Figure 5A
Row 4

Antibody incubation order: HBX → GAPDH → KI67 → OCT 3/4

LANE 1 2 3 4 5 6

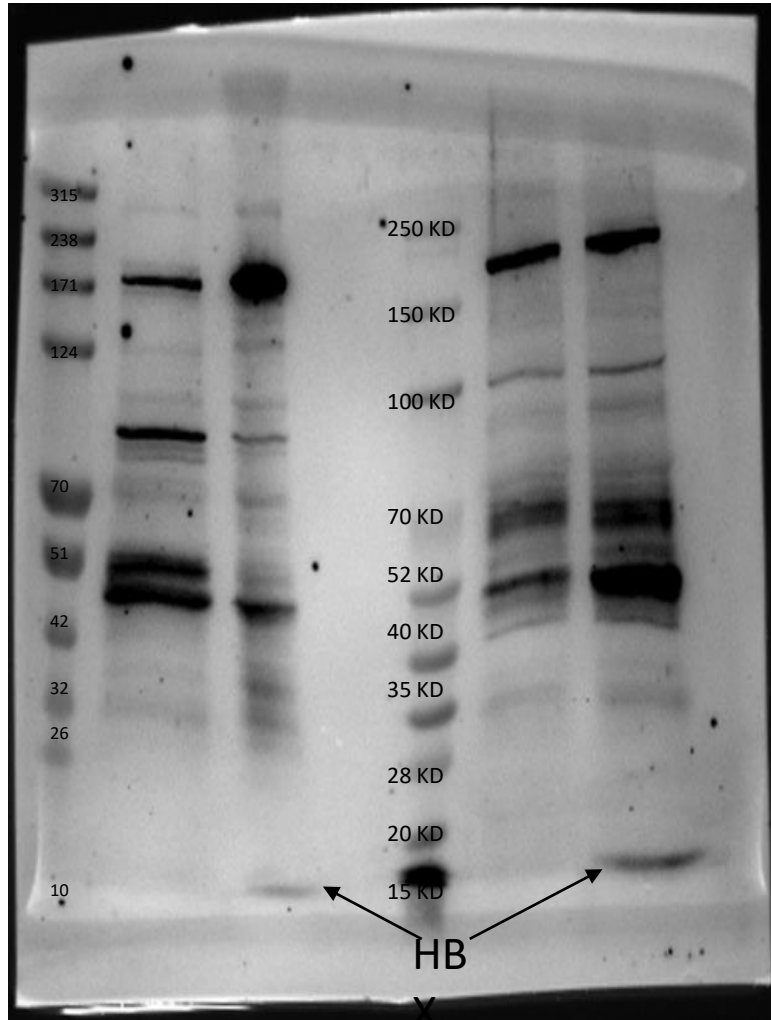

LANE

1

2

3

4

5

6

SAMPLES

MARKER

HUH7 Control

HUH7 HBX TF

MARKER

HEK Control

HEK pHBX TF

LANE

1

2

3

4

5

6



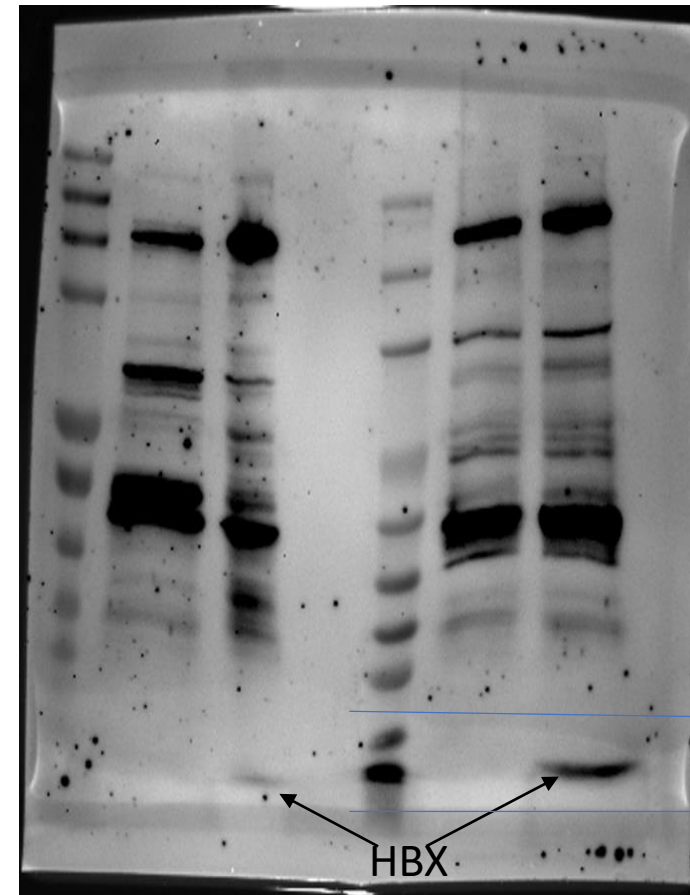


Figure-5A
Row 5

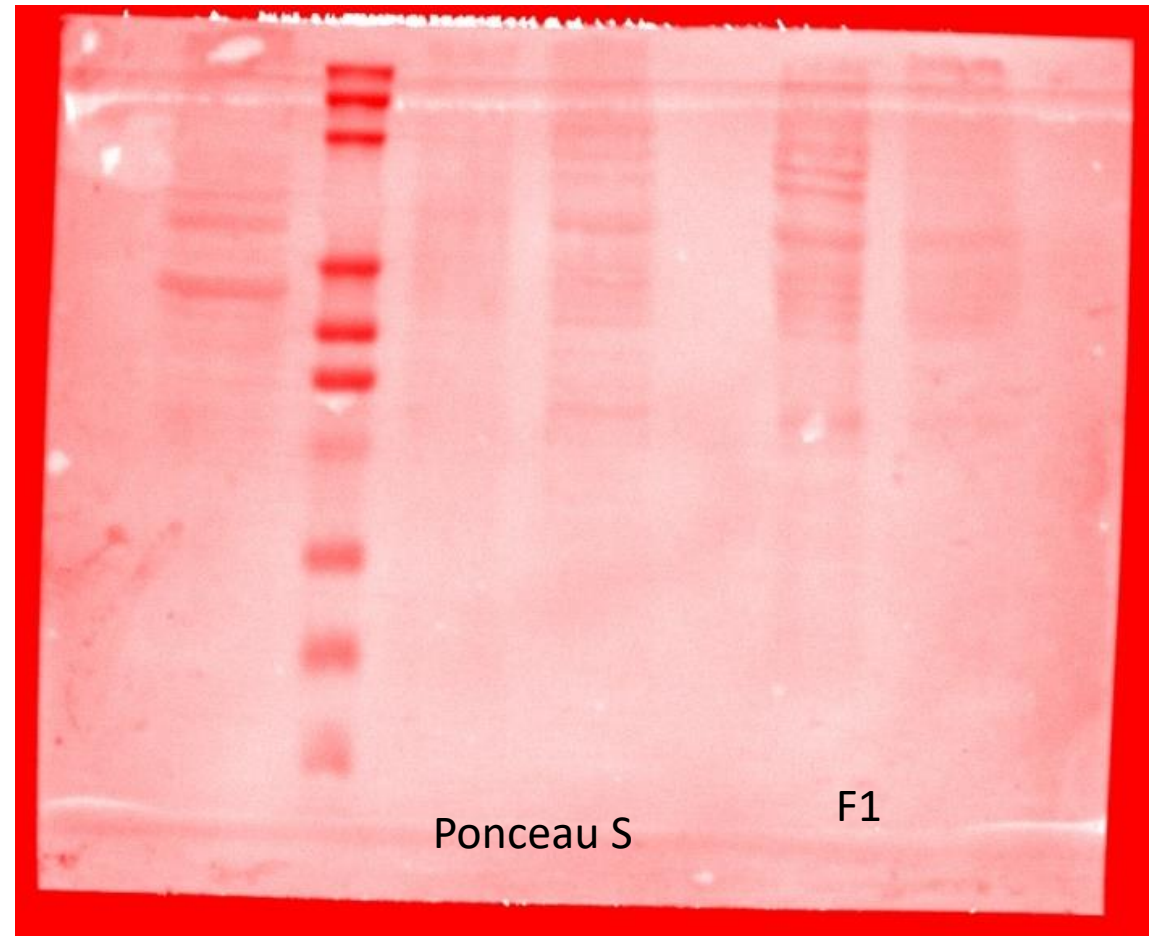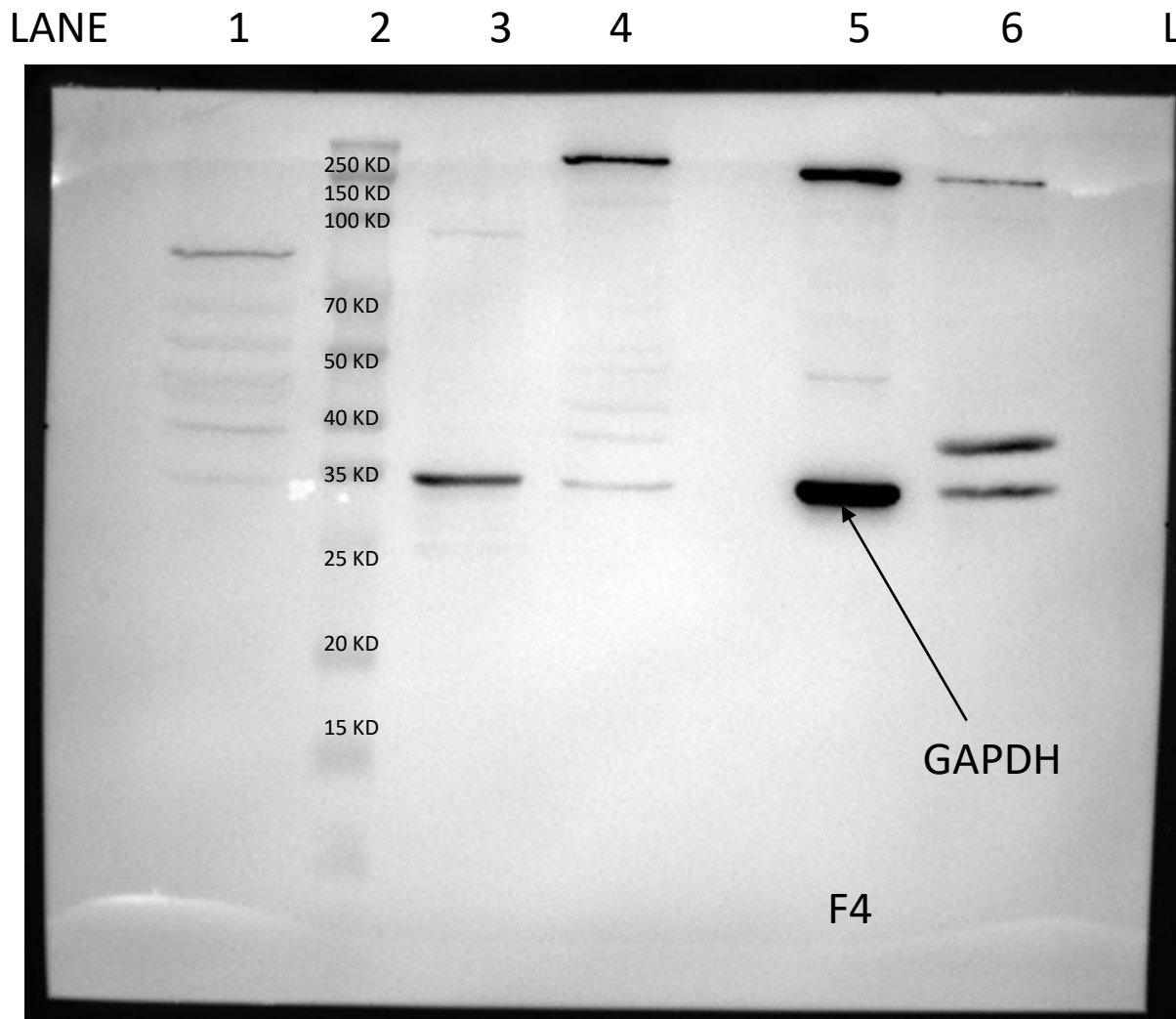

Antibody incubation order: OCT 3/4 → HBX → GAPDH

LANE	Samples
1	IPSC
2	MARKER
3	HEK CONTROL OLD
4	HEK HBX TF NEW
5	HEK CONTROL NEW
6	HEK HBX TF OLD

ANTIBODY	CAT NO	USED DILUTION	SIZE (aprox)
HBX	MA1-081	1:1000	17KD
OCT3/4	SC-5279	1:500	39KD
KI67	ITT2467	1:1000	367KD
GAPDH	60004-1-Ig	1:30000	36 KD
C-MYC	9402	1:1000	57 KD

WB Protocol:

1. 60 ug of total protein(BCA assay-read on Nanodrop) is loaded in each well of 12% SDS-PAGE.
2. Transfer at 90 V for 3 hrs.
3. Blocked with 5% BSA/5% Skimmed milk in TBST.
4. Overnight primary ab incubation.
5. Secondary ab incubation for 1-1.5 hrs.
6. Developed with WESTER ANTARE ECL kit.

OBSERVATIONS:

1. HBX is expressed in HEK Transfected sample.
2. KI67, OCT3/4 and c-MYC are overexpressed in HEK HBX Transfected sample.